



National Academy of Opticianry

## **Continuing Education Course**

Approved by the American Board of Opticianry

## **LENS POWERS**

Nominal Power-Focal Length  
Sagittal Depth  
Thick Lens / Back Vertex Power

### **National Academy of Opticianry**

8401 Corporate Drive #605

Landover, MD 20785

800-229-4828 phone

301-577-3880 fax

[www.nao.org](http://www.nao.org)

Copyright© 2016 by the National Academy of Opticianry. All rights reserved.  
No part of this text may be reproduced without permission in writing from the publisher.



# National Academy of Opticianry

## **PREFACE:**

This continuing education course was prepared under the auspices of the National Academy of Opticianry and is designed to be convenient, cost effective and practical for the Optician.

The skills and knowledge required to practice the profession of Opticianry will continue to change in the future as advances in technology are applied to the eye care specialty. Higher rates of obsolescence will result in an increased tempo of change as well as knowledge to meet these changes. The National Academy of Opticianry recognizes the need to provide a Continuing Education Program for all Opticians. This course has been developed as a part of the overall program to enable Opticians to develop and improve their technical knowledge and skills in their chosen profession.

The National Academy of Opticianry

## **INSTRUCTIONS:**

Read and study the material. After you feel that you understand the material thoroughly take the test following the instructions given at the beginning of the test. Upon completion of the test, mail the answer sheet to the National Academy of Opticianry, 8401 Corporate Drive, Suite 605, Landover, Maryland 20785 or fax it to 301-577-3880. Be sure you complete the evaluation form on the answer sheet. Please allow two weeks for the grading and a reply.

## **CREDITS:**

The American Board of Opticianry has approved this course for One (1) Continuing Education Credit toward certification renewal. To earn this credit, you must achieve a grade of 80% or higher on the test. The Academy will notify all test takers of their score and mail the credit certificate to those who pass. You must mail the appropriate section of the credit certificate to the ABO and/or your state licensing board to renew your certification/licensure. One portion is to be retained for your records.

## **AUTHORS:**

Diane F. Drake, LDO, ABOM, FCLSA, FNAO  
David F. Meldrum, LDO, ABOM, FNAO  
Randall L. Smith, M.S., ABOM, NCLEC, FNAO

## **INTENDED AUDIENCE:**

This course is intended for opticians of all levels.

## **COURSE DESCRIPTION:**

This course will explain lens power, nominal power and focal length of lenses. In addition, sagittal depth and thick lens design will be presented. Understanding these basic concepts helps the ECP to understand the final form of the patient's lens.

## **LEARNING OBJECTIVES:**

At the completion of this subsection, the student should be able to:

- Calculate nominal power
- Explain base curve and tool curves
- Calculate both focal length and dioptric power
- Be able to explain the relationship between focal length and lens power
- Discuss sagittal depth of a lens surface and how to calculate thickness of a lens
- Calculate the power of a thick lens

National Academy of Opticianry

8401 Corporate Drive #605

Landover, MD 20785

800-229-4828 phone

301-577-3880 fax

[www.nao.org](http://www.nao.org)

Copyright© 2015 by the National Academy of Opticianry. All rights reserved.  
No part of this text may be reproduced without permission in writing from the publisher.

**Lens Powers**  
**Nominal Power-Focal Length**  
**Sagittal Depth**  
**Thick Lens / Back Vertex Power**

Diane F. Drake, LDO, ABOM, FCLSA, FNAO  
David F. Meldrum, LDO, ABOM, FNAO  
Randall L. Smith, M.S., ABOM, NCLEC, FNAO


***Focal Length – The Definition***

Focal length is the point of convergence of light rays. In another definition, it is the distance between the back of the lens and the point focus.

The focal length of a lens is the distance at which the lens brings light to a focus.


***Positive Focal Length***

A plus lens creates a positive focal length, meaning that parallel light is brought to focus behind the lens. The image is real. The distance from the back surface of a plus lens to the point at which the light rays are brought to a point focus the positive focal length of that lens.

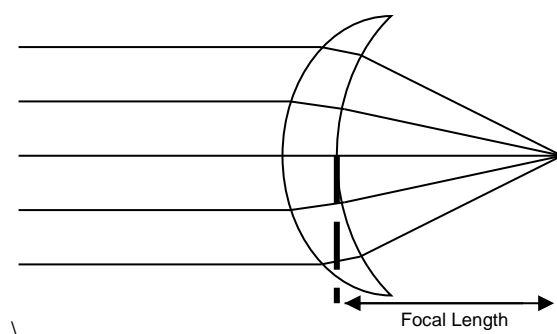

***Virtual Focal Length***

A minus lens creates a virtual focal length, meaning that the focal point is in front of the lens. The image is virtual. The distance from the back surface of a minus lens to the point at which the light rays would reach a common point if the diverged light rays were to be drawn back through the lens is called the virtual focal length for that lens.

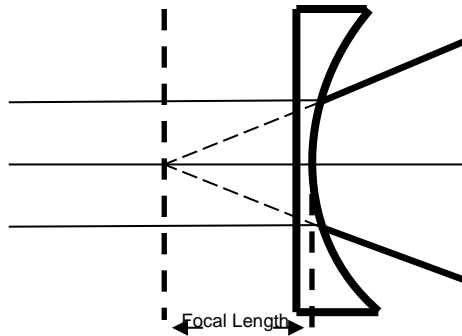

### *Dioptric Power of Lenses*

A lens diopter is defined as the lens power that will bring light to focus at a focal distance of 1 meter.

A diopter is a positive (or negative) optical unit of measurement, the reciprocal of the secondary focal length in meters, used to express the power of a lens. A plus (or minus) one diopter lens will converge (or diverge) parallel rays of light to a real (or virtual) focus at 1 meter.

The dioptric power of the lens has a direct, reciprocal relationship to the focal length. The stronger the dioptric power of the lens, the shorter the focal length, while the weaker the dioptric power, the longer the focal length. Other factors that influence dioptric power include thickness, radii of curvature and index of refraction which we will cover more in-depth later.

### *Focal Length*

The focal length in meters equals the reciprocal of the dioptric power of the lens. That means that if you know the dioptric power of the lens, to calculate the focal length, you simply divide the dioptric power into 1.

The formula is  $f = 1/D$

– In the formula:

- D = Dioptric power of lens
- f = focal length in meters

Remember to adjust formula based on mm, cm, dm or inches. In order to convert the formula for mm, cm, dm or inches you may convert by the metric system or simply use the formulae below. It's important for test purposes that you remember what formula you need to use.

- $f = 1/D = M = \text{meters}$
- $f = 10/D = \text{dm} = \text{decimeters}$
- $f = 100/D = \text{cm} = \text{centimeters}$
- $f = 1000/D = \text{mm} = \text{millimeters}$
- $f = 40/D$  (actually  $f = 39.37/D$ ) = inches

## *Lens Power*

The dioptric power of the lens equals the reciprocal of the focal length in meters. That means that if you know the focal length of the lens, to calculate the dioptric power of the lens, you simply divide the focal length into 1.

The formula is:  $D = 1/f$

– In the formula

- $D$  = Dioptric power of lens
- $f$  = focal length in meters

Again, remember to adjust formula based on mm, cm, dm or inches.

- $D = 1/f = M = \text{meters}$
- $D = 10/f = \text{dm} = \text{decimeters}$
- $D = 100/f = \text{cm} = \text{centimeters}$
- $D = 1000/f = \text{mm} = \text{millimeters}$
- $D = 40/f$  (actually  $D = 39.37/f$ ) = inches

## *Focal Length*

Example:

What is the focal length in cm of a +2.00 lens? Remember that we need to adjust formula for cm.

As in all optics, you will find there is more than one way to get to the same answer. Therefore, we will demonstrate two ways to find the answer. Choose the method that works best for you. Remember that we work in meters and the different formulas are simply variations on the meter.

The first method converts to centimeters AFTER finding the answer in meters. The second method converts to centimeters using the formula.

### *Method 1*

$f = 1/D$  using the formula, remember your answer will be in meters so you will have to convert to centimeters.

$$f = 1/D$$

$$f = \frac{1}{2} \text{ or } .5 \text{ meters}$$

.5 meters becomes  
5 decimeters becomes  
50 centimeters becomes  
500 millimeters

The answer in centimeters of a +2.00 D lens would be 50 cm.

### *Method 2*

$f = 1/D$  is the formula for meters, so the formula for centimeters is:

$$f = 100/D$$

therefore

$$f = 100/2$$

$$f = 50 \text{ cm}$$

### *Dioptric Power of a Lens*

Example:

What is the dioptric power of a lens with a virtual focal length (minus lens) of 400 millimeters?

Just as we demonstrated two methods for calculating the focal length, dioptric power can be found by using either of the two methods. Here are two examples.

The first method converts to diopters without changing the focal length to meters first. The second method converts the focal length to meters before using the formula.

### *Method 1*

$D = 1/f$  is the formula for meters so we have to adjust the formula for millimeters.

$$D = 1000/f$$

therefore

$$D = 1000/400$$



$$D = -2.50D$$

### Method 2

$D = 1/f$  is the formula for meters.

$$D = 1/f$$

Therefore, to convert 400 mm to meters


- .4 meters becomes
- 4 decimeters becomes
- 40 centimeters becomes
- 400 millimeters

$$D = 1/.4$$

$D = -2.50D$  (This is a minus power because it is a virtual focal length.)

Because we understand that every student comprehends in a different manner, we will present information in another format in addition to the one just discussed.

Here is a different approach to looking at focal lengths.

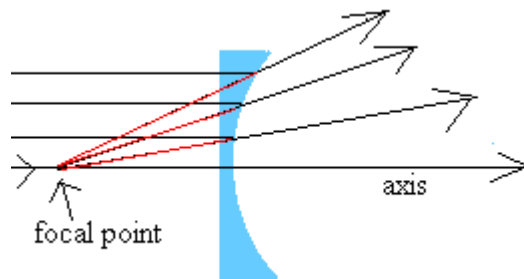

The lens pictured above has a virtual focal point (minus lens). The thickest part of the lens is on the edge and the thinnest part is in the center.

Light passing through a minus lens diverges, bends toward the base (lens edge). Light coming from a source closer than optical infinity (20 feet) is diverging. Light coming from a source at optical infinity or greater than optical infinity is considered parallel. If you look at the lens above, the light rays are diverging; therefore, the light doesn't come to a point focus behind the lens but appears to come to a point in front, thus a virtual focal point. The rays that emerge from this lens diverge from each other, instead of converging. We define divergence as negative vergence. This lens will add negative vergence to the rays entering it.


In this case, the rays emerge from the lens diverging *as if they had come from the focal point*, wherein the other lens they *actually crossed at the focal point*. In a lens that adds plus vergence

the focal point is *real* because the rays actually go through that point. In this case, the focal point is *virtual* because, although the rays act as if they came from the focal point, most of them never actually pass through the focal point.


You can start a fire by placing a piece of paper at the real focal point of the lens that adds plus vergence. This lens is a plus lens. You cannot start a fire with a virtual focal point, such as we have with the minus lens, which adds minus vergence.

You have probably used the word *virtual* before. "Virtual reality" or "virtual games" are all phrases that are used to describe something that seems very real, even though it is not. Well, a virtual image from a lens that adds negative vergence still appears to be very real; it just isn't!

When you look at an object, an image is formed on the light sensitive film (retina) in your eye, and your brain interprets the result as the object. The rays diverge from the object and are absorbed by the special light-sensitive nerve cells (rods & cones) in your eye. The processing that occurs in the nerve cells of the eye and brain interpret the object to be at the point that the rays diverged from.

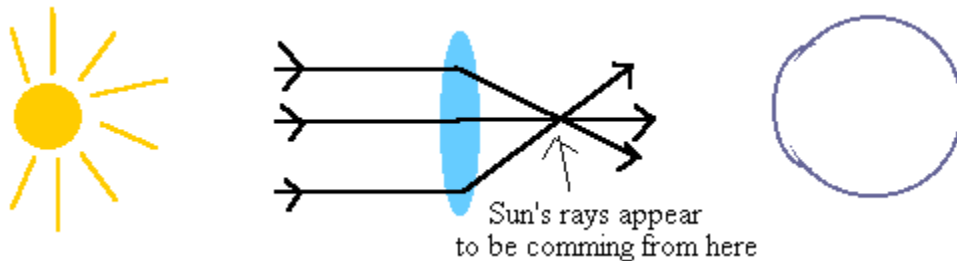

When a plus lens changes the vergence of the rays that come from a distance, the result is that the eye interprets the position of the object (in this case the sun) to be at the point that the rays are diverging from. (In these drawings the light is a long way away (beyond optical infinity), and the lens is several feet away from the eye. We are using the sun for illustration. We know that you know better than to look directly at the sun, through any lens or any other way.) **DO NOT TRY THIS AT HOME.**

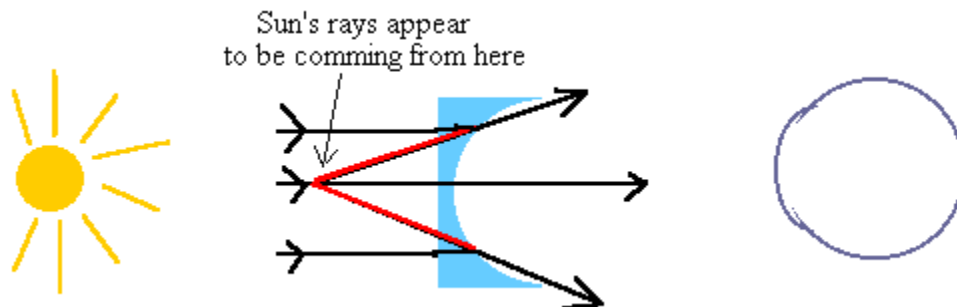

When a minus lens changes the vergence of the rays that come from the distance, the result is that the eye interprets the position of the object to be at the point where the rays *appear* to be diverging from, even though the rays never actually went through that point.

Well, what about the focal length?

With the plus lens, which added plus power, we had a real focal point, and the power of the lens was 1 divided by the focal length in meters. We will now say, by definition, that the focal length of a plus lens is positive, and the power in diopters is positive.

Then the lens that adds minus vergence is a minus lens, the focal length is negative, and the power is 1 divided by the negative focal length in meters, so the power of this lens is negative.

From this moment on, whenever you specify a power for a lens, it must have a + or a - sign. If it does not have a sign it will be wrong; whether the power is right or not!

Suppose you have a lens with a virtual focal length of 1 meter. What is the lens power?

We are going to use  $f$  to represent focal length and  $D$  to represent diopters. A virtual focal length is negative, so the focal length is -1 m.  $D = 1 \div f = 1 \div (-1)$ , so  $D = -1$  diopters. This is another way to write the same formula that we used previously.

$D = 1/f$  is the same formula as  $D = 1 \div f$

If the virtual focal length is 0.5 m, then the power is  $D = 1 \div f = 1 \div (-0.5)$ , so  $D = -2.00$  diopters.

A minus lens has a focal length of 40 cm. What is its power?

A minus lens has a negative focal length, so  $f = -40$  cm.

The focal length must be in meters, so  $f = -40$  cm = -0.40 m.

$D = 1 \div f = 1 \div (-0.4)$  m

$D = -2.50D$

Right?

What is the focal length of a -6.25D lens in mm?

$D = -6.25$

$f = 1 \div D = 1 \div (-6.25)$

$f = -0.16$  m = -160 mm. (That will be a virtual focal point in front of the lens.)

Make sure that you keep track of the minus sign and the decimal point!

We have discussed that 40 inches are almost equal to one meter. So, if you are given the focal length in inches, or if you are given the power and asked for the focal length in inches, you can use the equation  $D = 40 \div f$  or  $f = 40 \div D$ .

In this case, the negative lens (minus lens) still has a negative focal length, and the positive lens (plus lens) still has a positive focal length.

### ***Nominal Power of Lenses***

By definition, the word nominal means in name only. That means that it is not exact, or actually, it is in name only. A real-life example of nominal would be that if you go to the lumber yard to buy 2x4's, what will you get? Lumber with dimensions of 2" by 4"? You get lumber with dimensions 1.75" x 3.50". 2x4 is the nominal dimension of the lumber, but the lumber is not, in fact, that size. So when we are talking about the power of a lens, the nominal will be an exact power but will be a very close approximation.

The definition of the nominal power of a lens is the sum of the front and the back surface powers.

The Nominal Power formula is:

$$D_e = D_1 + D_2$$

Where:

$D_e$  = effective, vertex or lensometer power

$D_1$  = true power of front curve

$D_2$  = true power of back curve

As has been discussed previously, you may encounter varying versions of a formula. Another form of the Nominal Power formula is:

$$D_N = D_1 + D_2$$

Where:

$D_N$  = effective, vertex or lensometer power

$D_1$  = true power of front curve

$D_2$  = true power of back curve

Here we have the nominal power formula for our lenses. What does this mean? Does it mean that the lens will not really be the power that we get from this formula?

That is exactly what it means!

All of the following are the variables that determine the actual power of a lens:

- The curvature of the front surface
- The curvature of the back surface
- The material the lens is made from
- The material on each side of the lens (may be the same on each side, or may be two different materials)

- The thickness of the lens
- The wavelength of the light going through the lens
- Which surface the ray of light enters first

For our purposes, knowing the nominal power of a lens will assist us in frame and lens material selection. Items 3 – 7 are higher order lens parameters that will affect the actual power of the lens.

The *Nominal Power* formula uses just the surface curves of the lens. Some people call this the thin lens power formula; because it assumes that the lens has no thickness, so the thickness will not affect the power of the lens. This nominal power formula is fine for lenses with a power that is less than about +4.00 D; above that the assumption that the lens thickness is negligible is no longer valid.

We used curvature and lens material index to find the power of each surface. Now we will find the lens power from the surface powers.

The nominal lens formula says that you can add together the powers of the front and back surfaces of the lens, and you will have the nominal power of the lens.

When we defined surfaces, we defined a concave (back) surface as having minus power and a convex (front) surface as having plus power.

Let's look at some examples.

$$D_e = D_1 + D_2$$

$D_1$  = Front (base) curve (power of the front surface)

$D_2$  = Back (ocular or tool) curve/s

1. Find the nominal power ( $D_e$ ) if:
  - $D_1 = +10.00$  and  $D_2 = -6.00$  - see sample below
  - $D_1 = +2.00$  and  $D_2 = -4.50$
  - $D_1 = +8.00$  and  $D_2 = -5.75$
  - $D_1 = -1.25$  and  $D_2 = -3.75$
  - $D_1 = -5.00$  and  $D_2 = -6.00$

Example:

$$D_e = D_1 + D_2$$

$$D_e = (+10.00) + (-6.00)$$

$$D_e = +4.00$$

Same formula format:

$$D_e = D_1 + D_2$$

$D_1$  = Front (base) curve (power of the front surface)

$D_2$ = Back (ocular or tool) curve/s

2. Find the back surface ( $D_2$ ) if:

- $D_e = +10.00$  and  $D_1 = +6.00$  – see sample below
- $D_e = +2.00$  and  $D_1 = +8.50$
- $D_e = +8.00$  and  $D_1 = +5.75$
- $D_e = -5.00$  and  $D_1 = -1.00$
- $D_e = +6.50$  and  $D_1 = +10.00$

Example:

$$\begin{aligned}D_e &= D_1 + D_2 \\+10.00 &= +6.00 + D_2 \\+4.00 &= D_2\end{aligned}$$

Now let's give some examples using the different form for the formula:

$$\begin{aligned}D_N &= D_1 + D_2 \\D_1 &= \text{Front (base) curve (power of the front surface)} \\D_2 &= \text{Back (ocular or tool) curve/s}\end{aligned}$$

3. Find the front surface ( $D_1$ ) if:

- $D_N = +10.00$  and  $D_2 = -2.00$
- $D_N = +2.00$  and  $D_2 = -8.50$
- $D_N = +8.00$  and  $D_2 = -2.75$
- $D_N = -1.25$  and  $D_2 = -7.75$
- $D_N = -5.00$  and  $D_2 = -3.00$
- $D_N = +6.50$  and  $D_2 = -10.00$

$$\begin{aligned}D_N &= D_1 + D_2 \\+10.00 &= D_1 + (-2.00) \\+12.00 &= D_1\end{aligned}$$


4. I need a lens with power +3.50D. If I use a lens blank with a front surface of +8.00D, what back surface do I need to make?
5. I have a lens with a front surface of +9.25D and a back surface of -2.25D. What is the power of the lens?
6. I want my lens to have a power of -3.00D. What back surface do I need if the front surface is +8.00D?

You may encounter the terminology “bent design lenses”. This simply means that you have one surface that is plus/convex and the other surface is minus/concave. In a conventional cylindrical lens, you will have two curves on the back with one on the front, while a conventional spherical

lens only has one curve on the back one curve on the front. More technologically advanced surfaces and designs will be covered later in the program.

### ***Sagittal Depth/Sagittal Value***

The definition: Sagittal depth is the straight line distance between the back surface of a lens at its vault and the chord diameter of the lens. The height (or depth of a convex or concave lens surface is referred to as the sagitta or more oftentimes simply SAG of the surface. Simply put, it is the perpendicular distance from the vertex of the curve to some plane cutting through the curve. It is also referred to as the vertex depth of the lens surface. The diameter of a lens is calculated by the distance (or chord) from one side of the curve to the other across the plane that cuts through the curve. When we use a lens clock (lens measure) in effect we are measuring the SAG of the curve. We translate that to the power of the curve. For a spherical surface, the sagitta of a lens surface is defined by the radius of curvature of the surface at a given diameter (that is, the size of the lens blank).

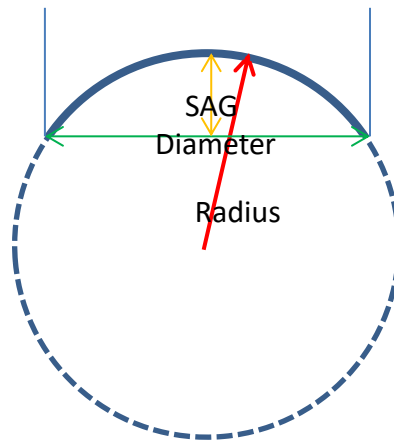


In addition, the sagitta of a *convex* curve is equivalent to the sagitta of a *concave* curve when the radii of curvature of both surfaces are equal. So what is the significance of sagittal depth? We can calculate the center thickness of a convex lens or the edge thickness of concave lens with the use of the sagittal depth formula. The formula is:


- Thickness =  $\frac{\text{Radius}^2 \times \text{Power}}{2000 (N - 1)}$

- Where:
- N = index of refraction of lens material used
- Thickness is in mm

The illustration below demonstrates how the surface curvature and the diameter relates to the sagittal depth/height of the lens. In determining actual thickness, keep in mind that the formula must be calculated and then added to the edge thickness of a plus lens or the center thickness of a minus lens. *The formula is based upon 0 CT or edge thickness.*

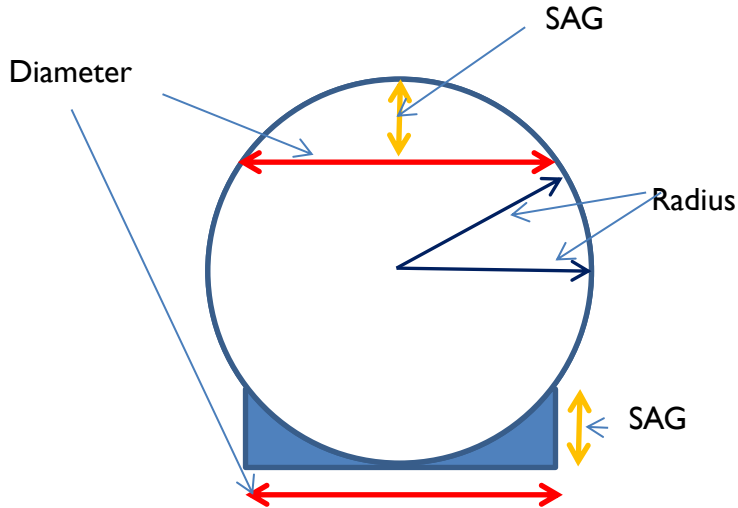

Another illustration that may help you visualize SAG is below.

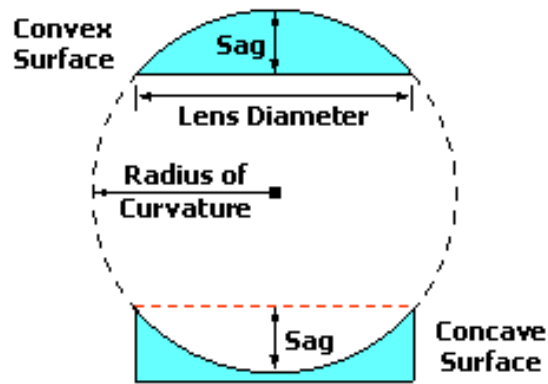


**Geometry of a Lens Surface**



The following illustrations may help you to visualize the relationship.


### Sagittal Depth Relationship

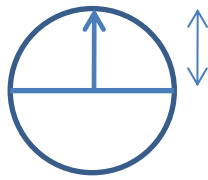


One radius




Shorter radius

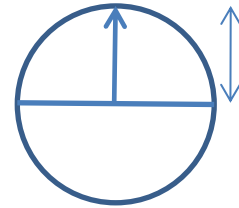


Larger radius

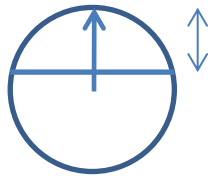

One radius

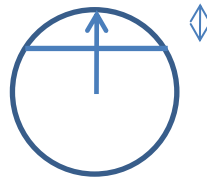




Same radius  
Shorter diameter  
Shorter sagittal depth

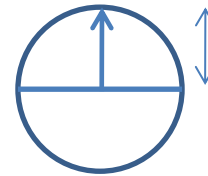

Same radius  
Longer diameter  
Longer sagittal depth


### Same Radius - Different Diameters/Different Sagittal Depths

If the radius remains the same and the diameter changes, the thickness will change, and the sagittal depth will change. For example, a shorter radius will create a steeper surface curvature. A longer radius will create a flatter surface curvature. If the radius remains the same, but the diameter is smaller, the sagittal depth will be shorter. If the radius remains the same but the diameter is larger, the sagittal depth will be longer.

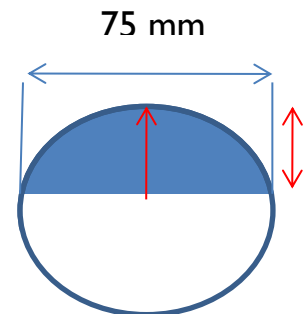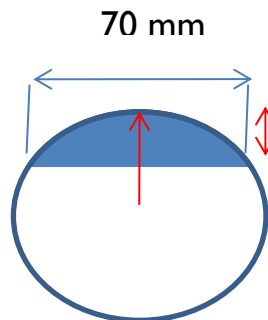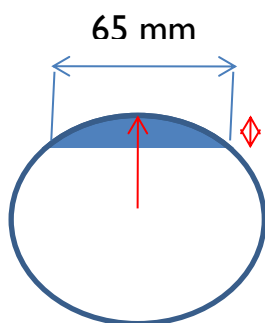


### Same Sagittal Depth - Different Radius/Different Diameter

If the sagittal depth remains the same, but the diameter increases, the radius will also increase. As the diameter and the radius increase, the surface becomes flatter.

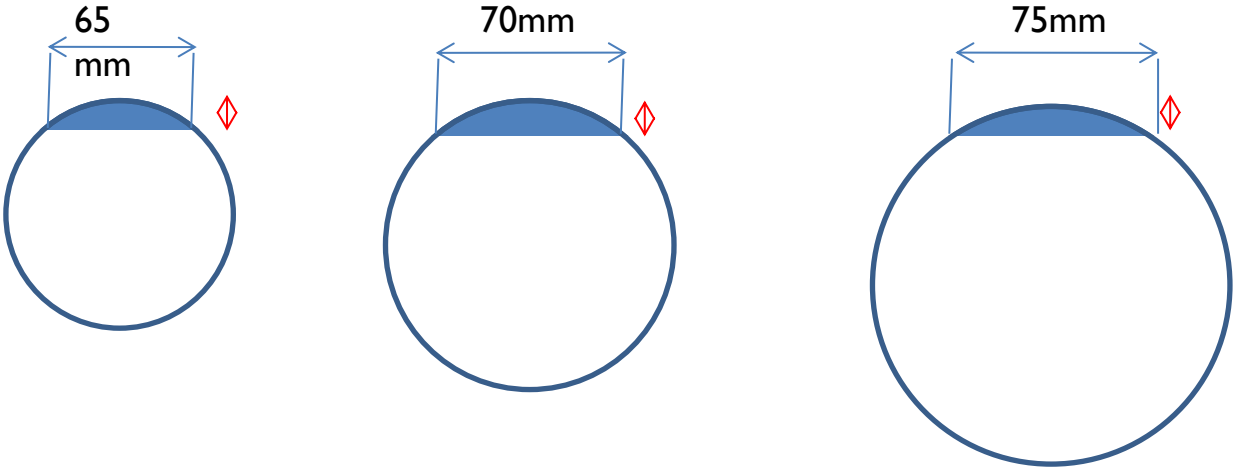


### Same Diameter - Different Radius/Different Sagittal Depths

If the diameter remains the same, but the radius changes, the sagittal depth will also change. For example, if the diameter remains constant, but the radius increases, the sagittal depth will also be smaller, creating a flatter surface.

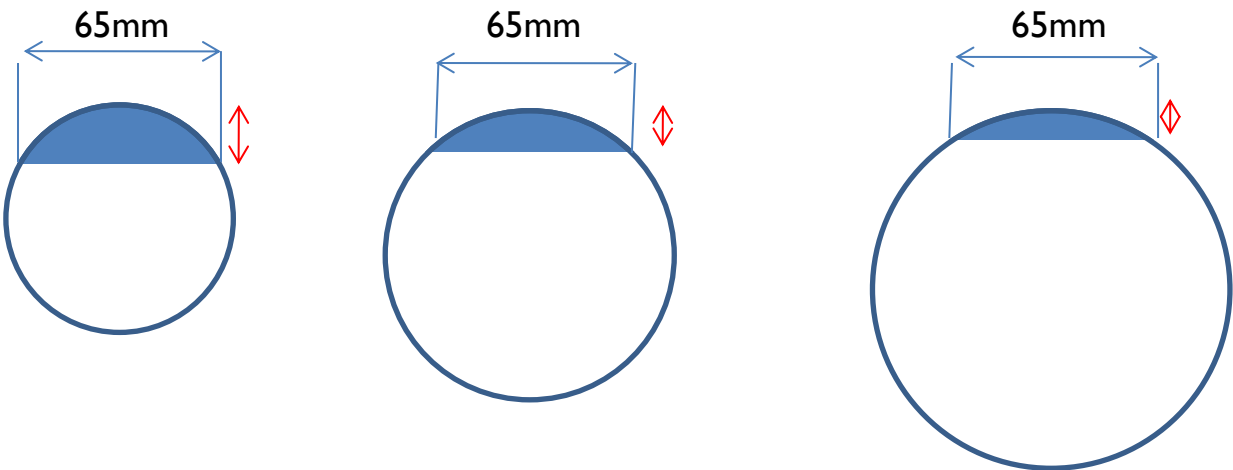

The illustration below demonstrates how a lens clock/lens measure measures sagittal depth.

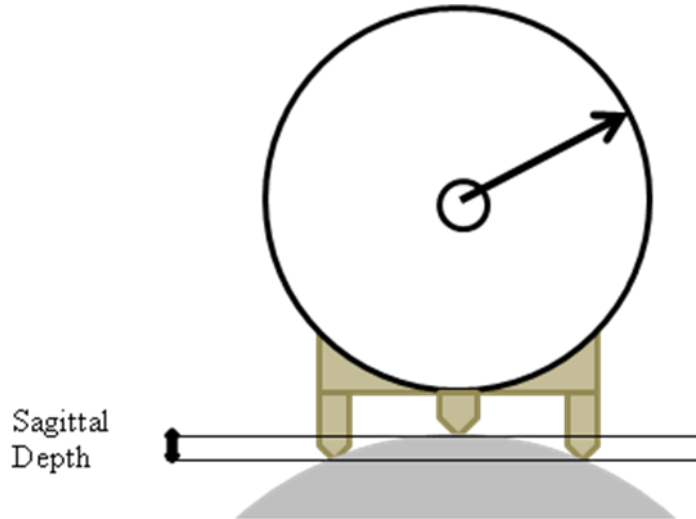

### ***Thick Lens / Back Vertex Power***

The back surface or ocular surface of a lens sits closest to the eye. Therefore it is important that we be able to calculate the power of the lens based upon the position in front of the eye. In neutralizing the power of a lens on a lensometer, we position the back surface of the lens against the lens stop. That gives us the back vertex power. In thinner lenses, there is very little difference between front surface power and back surface power as read in a lensometer. However, as the thickness of the lens increases, there can be a significant difference between the two readings.

Remember that the thin lens/nominal lens formula goes on the theory that there is NO thickness, so thickness has no influence on the power.


Vergence: The definition - The amount of divergence (negative) or convergence (positive) of a pencil of light rays entering or leaving a lens.

Based upon the principles of optics, as light strikes an optical medium at an oblique angle, the ray is bent toward the normal as it is passing through; which means it is bending (either divergent or convergent), within the lens.

For thick lens designs, we must also consider the equivalent thickness ( $t/n$ ), which is the vergence of the light passing through the thickness of the lens. Because of the steep front curve and the thickness of the center, the thick lens will add convergence to the wave front as it passes through the lens. To find the back vertex power (BVP) in addition to the simple combination of the powers of the front and back surfaces, we must also include the factors of the thickness of the lens in meters and the index of refraction of the material.

For example, using the nominal power/thin lens formula we would calculate the power and it would appear as the illustration directly below.

#### Thin Lens Design

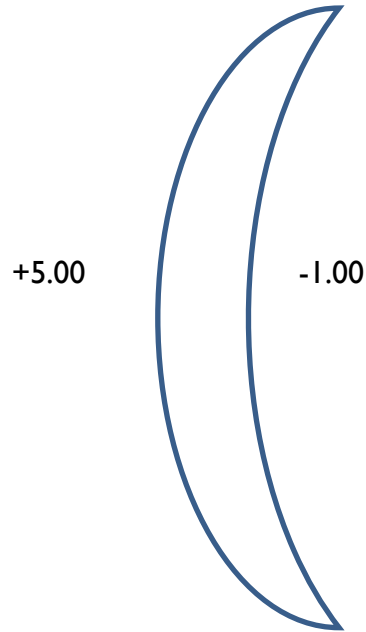

Approximate power +4.00

However, when the power of a lens is over +/- 4.00D or if a steeper base curve is used, we need to consider the thickness as well, due to vergence within the lens.

Therefore, for the thick lens formula, we must include the equivalent thickness ( $t/n$ ) to the formula.

The thick lens/back vertex formula is:

- $D_e = D_1 + D_2 + (t/n \times D_1^2)$


Where:

- $D_e$  = effective, vertex or lensometer power
- $D_1$  = true power of front curve
- $D_2$  = true power of back curve
- $t$  = lens thickness in meters
- $n$  = index of refraction

The illustration below demonstrates the thick lens design. If we simply used the thin lens formula, the power would appear to be +4.00D, however, due to the high surface curvature, the

thickness is increased, which not only increases the magnification effect of the lens but the power as well. Therefore, the approximate power using the thick lens formula would be +4.39D.

### Thick Lens Design

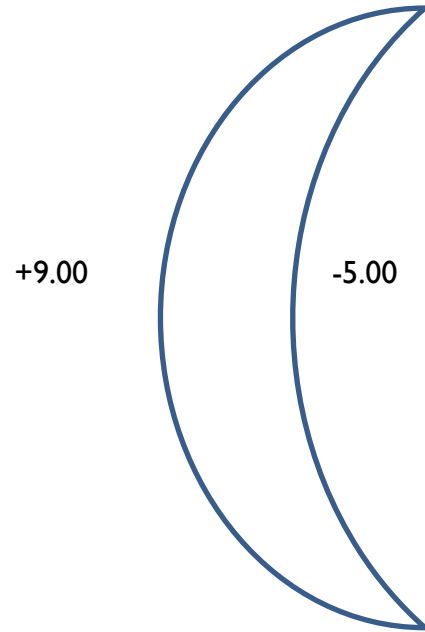

$t = 7.0\text{mm}$   
 $n = 1.50$   
Approximate power +4.39

National Academy of Opticianry

8401 Corporate Drive #605

Landover, MD 20785

800-229-4828 phone

301-577-3880 fax

[www.nao.org](http://www.nao.org)

Copyright© 2015 by the National Academy of Opticianry. All rights reserved.  
No part of this text may be reproduced without permission in writing from the publisher.